
Camry

— Premium —



CR 8185

- | | |
|--------------------------------------|-------------------------------------|
| (GB) user manual - 3 | (DE) bedienungsanweisung - 4 |
| (FR) mode d'emploi - 7 | (ES) manual de uso - 9 |
| (PT) manual de serviço - 11 | (LT) naudojimo instrukcija - 13 |
| (LV) lietošanas instrukcija - 15 | (EST) kasutusjuhend - 17 |
| (HU) felhasználói kézikönyv - 22 | (BS) upute za rad - 20 |
| (RO) Instrucțiunea de deservire - 18 | (CZ) návod k obsluze - 28 |
| (RUS) инструкция обслуживания - 30 | (GR) οδηγίες χρήσεως - 24 |
| (MK) упатство за корисникот - 26 | (NL) handleiding - 33 |
| (SL) navodila za uporabo - 35 | (FI) manual ng pagtuturo - 38 |
| (PL) instrukcja obsługi - 62 | (IT) istruzioni operative - 44 |
| (HR) upute za uporabu - 36 | (SV) instruktionsbok - 40 |
| (DK) brugsanvisning - 48 | (UA) інструкція з експлуатації - 50 |
| (SR) Корисничко упутство - 46 | (SK) Používateľská príručka - 42 |
| (AR) دليل التعليمات - 52 | (BG) Инструкция за употреба - 54 |
| (AZ) istifadəçi kitabçası - 56 | (ALB) Udhëzime për përdorim - 58 |

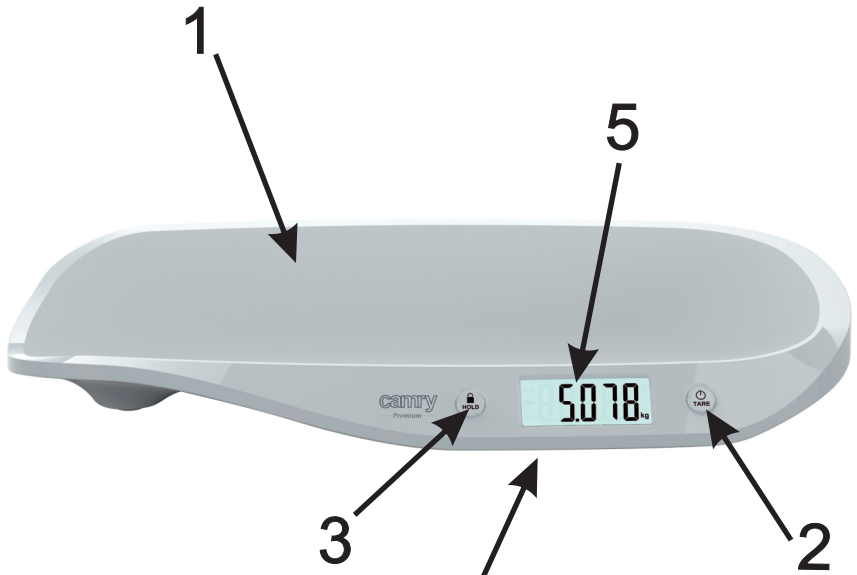4

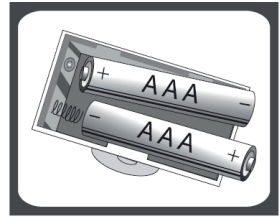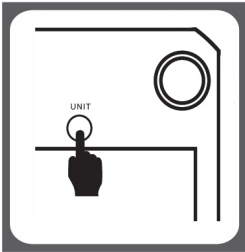a

Lo

b

Eeeee

SAFETY CONDITIONS. IMPORTANT SAFETY INSTRUCTIONS PLEASE READ CAREFULLY AND KEEP FOR FUTURE REFERENCE

The device is for home use only.

1. Before using the product, read the instructions below carefully and always follow them. The manufacturer is not responsible for any damage resulting from improper use of the product. Do not use for any other purposes inconsistent with its intended use.
2. The product is intended for indoor use only.
3. Be careful when using the product when children are nearby. Do not let children play with the product. Do not allow children or people who are unfamiliar with the device to use it unsupervised.
4. **WARNING:** This appliance can be used by children over 8 years of age and persons with reduced physical, sensory or mental capabilities or lack of experience or knowledge of the appliance only if they are supervised by a person responsible for their safety or have been instructed on the safe use of use of the device and are aware of the risks associated with its operation. Children should not play with the device. Cleaning and maintenance of the device should not be performed by children unless they are over 8 years of age and these activities are performed under supervision.
5. Never immerse the entire device in water. Never expose the product to atmospheric conditions such as direct sunlight, rain, etc. Never use the product in humid conditions. Never use the device near water, e.g. near a shower, bathtub or swimming pool.
6. Never use the product if it has been dropped or otherwise damaged, or if it is not working properly. Do not attempt to repair a damaged product yourself. Always contact an authorized service center for repair if your device is damaged. All repairs may only be performed by authorized service technicians. Repairs performed incorrectly may pose a danger to the user.
7. Never place the product on or near hot or warm surfaces or kitchen appliances such as an electric oven or gas burner.
8. Never use the product near flammable materials.
9. Do not touch the device with wet hands.
10. The batteries powering the device may leak if they are discharged or have not been used for a long time. To protect your device and your health, replace it regularly and avoid skin contact with leaking battery electrolyte.
11. If the scale will not be used for a long time, remove the battery.
12. Unless the device has an automatic shut-off function, it should be turned off manually after each use.
13. Use the device on a stable and flat surface.
14. Protect the device from the influence of strong electromagnetic fields.
15. Make sure the surface on which the device is used is not wet or slippery.
16. Be careful not to drop heavy objects on the scale's display or glass.
17. **DO NOT** hit, shake or drop the scale.
18. During the entire measurement, you must stand still on the scale (applies to scales for

measuring body weight)

DEVICE DESCRIPTION

1. Scale surface	2. ON/OFF and TARE button	3. Hold button
4. UNIT button	5. LCD display	6. Battery compartment

BEFORE FIRST USE

Insert the correct batteries into the compartment (See "Replacing the battery")

REPLACING THE BATTERY

Remove the lid of the battery compartment at the bottom of the scale. Remove the old battery. Insert the new battery of the right type. Pay special attention to the polarity (+) plus and (-). Close lid of the battery compartment.

SELECTING UNITS OF MEASUREMENT

Press and hold the ON/OFF button (2) to turn on the scale. Press the UNIT button (4) to change the units of measure. UNIT button (4) is placed at the bottom of the scale. The metric system will change to kg, lb or oz. The unit of measurement will be visible on the display (5).

WEIGHING A CHILD

Press and hold the ON/OFF button (2) to turn on the scale. Wait until the scale displays the 0.0 value.

Carefully place the child on the surface of the scale (1). Be careful that the baby doesn't fall. Pay special attention to the baby's head.

Wait a moment for the result to stabilize.

NOTE: Do not touch the baby while the scale takes the measurement, as this can affect the final result.

As soon as the baby's weight is determined, the word "HOLD" on the display (5) will stop flashing.

Display will show result for 20 seconds.

HOLD and Auto-HOLD function

Scale has HOLD and Auto-HOLD function. The Auto-HOLD function will show the baby's weight on the display (5) as soon as baby holds still enough for an exact measurement to be made.

If the child is moving all the time, and the result cannot be stabilized, press HOLD button (3).

Within 3 seconds the word "HOLD" will stop flashing on the display (5) and the result will stabilize.

Display will show the result for 20 seconds. Then the scale will switch off automatically.

TARE function

The TARE function allows for weighing the child on a towel, blanket or diaper.

For example, place a blanket on the scale and press TARE button (2). Once you press the TARE button (2) "TARE" sign will be displayed on the screen. Next, place the child on the scale.

CLEANING AND MAINTENANCE

– Gently wipe the scale with a dry or slightly damp cloth.

MESSAGES

EEEEEE – excessive load

LO – low battery

TECHNICAL DATA

Power: 2x 1,5V (AAA)

Maximum load: 20kg / 44lb / 705oz

Minimum load: 0.05kg / 0.01lb / 0.18oz

Please transfer cardboard packaging to waste paper. Pour the polyethylene bags (PE) into the plastic container. Worn device should be returned to the appropriate point storage because it is dangerous in the device ingredients may pose a threat to the environment. The electrical device should be handed over to limit its reuse and use. If in the device there are batteries you have to take out and give to the point stored separately.

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